

# specifications

The medical gas outlet shall be an Amico Alert-1 series.

Outlets shall be manufactured with a 6-1/2" [165 mm] length type "K" 1/2" [12.7 mm] outside diameter (3/8" nominal) size copper inlet pipe stub, which is silver brazed to the outlet body. The body shall be of 1-5/16" [33 mm] diameter, one piece brass construction. For positive pressure gas services, the outlet shall be equipped with a primary and secondary check valve and the secondary check valve shall be rated at a maximum of 200 psi [1,379 kPa] in the event the primary check valve is removed for maintenance. Outlet bodies shall be gas specific by indexing each gas service to a gas specific dual pin indexing arrangement on the respective identification module.

A large color coded front plate shall be used for ease of gas identification and aesthetic appeal.

A one piece chromed fascia plate shall frame the outlet. With the back rough-in mounted, the outlet shall adjust from 3/8" [10 mm] to 1" [25 mm] variation in wall thickness.

The outlets shall be of modular design and include a gas specific 16 ga. [1.6 mm] steel mounting plate designed to permit on-site ganging of multiple outlets, in any order, on 5" [127 mm] spacing.

The latch-valve assembly shall be Puritan-Bennett Quick Disconnect compatible and accept only corresponding Puritan-Bennett gas specific adapters.

All M.R.I. outlets are manufactured from nonferrous materials.

All outlets shall be cleaned and degreased for medical gas service, factory assembled and tested.

The outlets are UL listed.

Amico products comply with NFPA-99.



# features

- Flow rate exceeds the requirements of NFPA and CSA
- Accepts only Puritan-Bennett gas specific adapters
- 100% Hydrostatic testing of each outlet and cleaned for medical gas service
- Indexed to eliminate interchangeability of gas services
- Up to 1" [25 mm] wall thickness adjustment
- 360° Swivel Inlet Pipe for easy installation
- Universal Rough-In Assembly to accept Quick Disconnects or DISS Front Adapters (Can be interchanged at any time)
- Outlets are also available for M.R.I. use

**project**

---




---



---

# wall outlet quick disconnect - puritan-bennet compatible

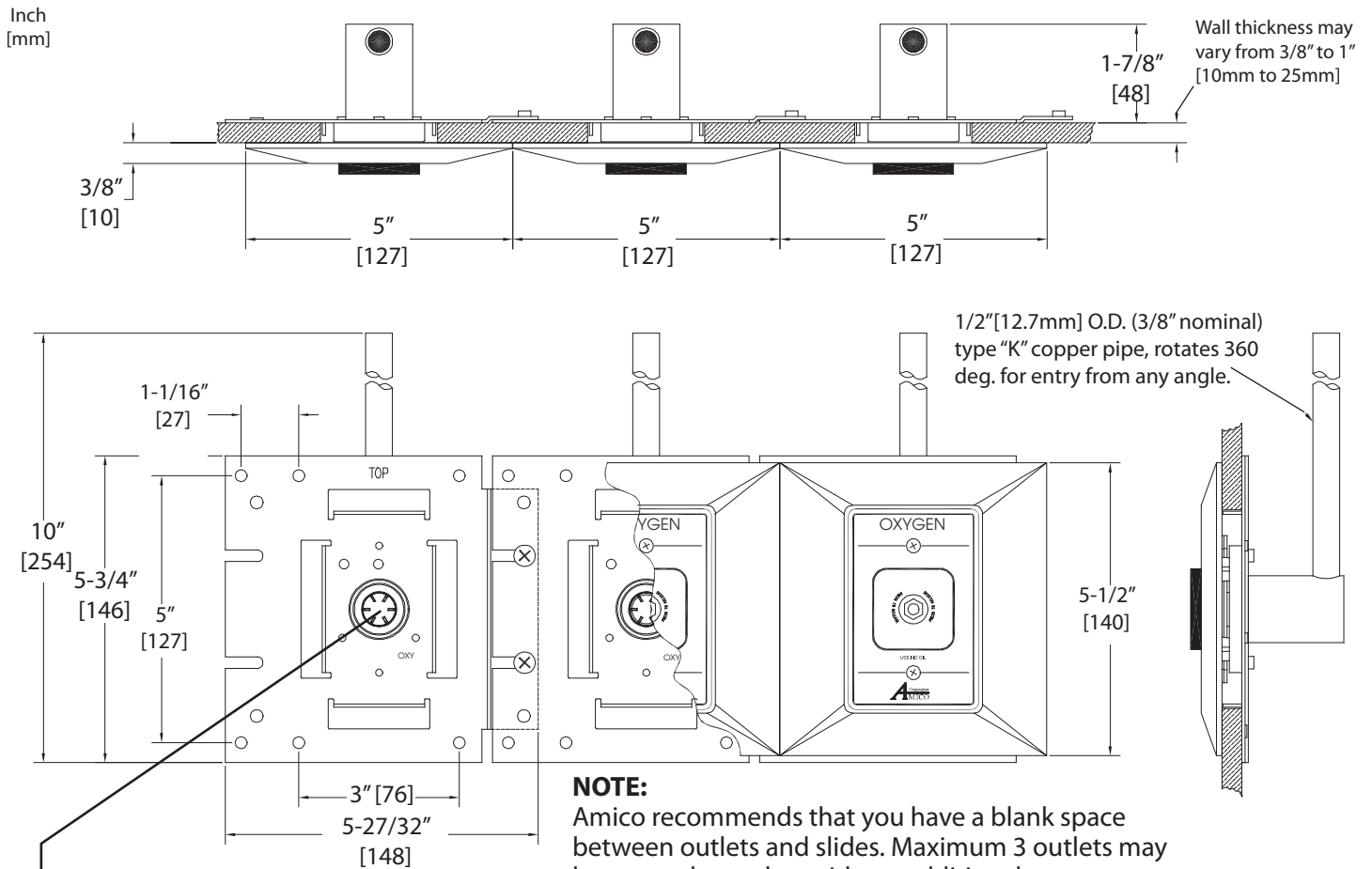

<sup>1</sup>NOTE: For CSA/ISO, suction inlets are not supplied with a secondary check valve. A pressure plug is provided for testing purposes, rated at 150 psi [1,034 kPa].

**NOTE:**

Amico recommends that you have a blank space between outlets and slides. Maximum 3 outlets may be ganged together without additional support.

## model numbers

### O-PBWAL-L-XXX

The L defines the Language:

- U = English (NFPA) Oxy, Air, Vac and WAG only
- E = English (CSA/ISO)
- F = French (CSA/ISO)
- S = Spanish (NFPA)

Example: Oxygen NFPA = **O-PBWAL-U-OXY**

## M.R.I. model numbers

### O-PBWALM-L-XXX

Example: English NFPA (M.R.I.) - Puritan-Bennett, Oxygen = **O-PBWALM-U-OXY**

The XXX defines the Gas:

- Oxygen = OXY
- Medical Air = AIR
- Nitrous Oxide = N2O
- MedVac<sup>1</sup> = VAC
- Waste Anesthetic Gas Disposal (NFPA) = WAG
- Anesthetic Gas Scavenging System<sup>1</sup> (CSA) = AGS

<sup>1</sup>Please see note above

Oxygen Systems Plus, Inc.  
1-800-480-7664

represented by: